



-  Whitening
-  Skin Rejuvenation
-  Pore Reduction
-  Acne Scar
-  Wrinkle Reduction

Technical Specifications

| | |
|----------------------------|------------------------------|
| Model name | VIVACE |
| Power Input | 100/230VAC 50/60Hz |
| Electric power consumption | 40VA |
| Dimensions | 380(W) X 1140(H) X 340(L)mm |
| Weight | 27g |
| output | |
| Frequency | 1Mhz / 2Mhz (High Frequency) |
| output power control | PWM |
| interface | 10.4 inch LCD Touch screen |

VIVACE

TRUE FRACTIONAL MICRONEEDLE RF

Dermal Volumetric Rejuvenation(DVR) System:
Fractional Micro-needle RF





Wrinkle Reduction



Acne Scar



Pore Reduction



Skin Rejuvenation



Whitening

SMART & SIMPLE VIVACE

Distinction of UIVACE

- Micro-needle electrodes deliver RF energy directly into the dermis
Thermal injury induces the production of new collagen and elastin
- Minimally invasive dermal volumetric rejuvenation system
less pain & down time and quick results
- Special golden & partly insulated Micro-needle target the dermis
without epidermal injury
- Precise 1MHz / 2 MHz RF makes dramatic effect
- Softly inserting Micro memory motor needling reduce pain & adverse effect
- Well organized controls for the various skin conditions
- Save Program for easy use
- Red & Blue Light emitting diode(LED) lights help skin activity and acne
improvement
- User friendly and safe Single use cartridge

Accessory introduction

- 36 gold coated micro-needles
- Step-Motor System
- User Friendly ergonomic & lighter applicator
- Touch screen
- Save programs
- Hand & Foot Dual Switch
- Level 1-9 / Depth 0.5-3.5mm





- ✓ Quick Result!
- ✓ Less Pain!
- ✓ Less Downtime!



Science of UIVACE

1. CIT(Collagen Induction Therapy) by Micro-needling

When it is applied to the skin, under local or topical anesthesia, sterile Micro-needle are used to create many microscopic channels deep into the dermis of the skin, which stimulate your own body to produce new collagen. These channels also improve the penetration of vitamins A and C creams. Vitamins A and C stimulate skin renewal, thereby making the skin appear fresher and younger.

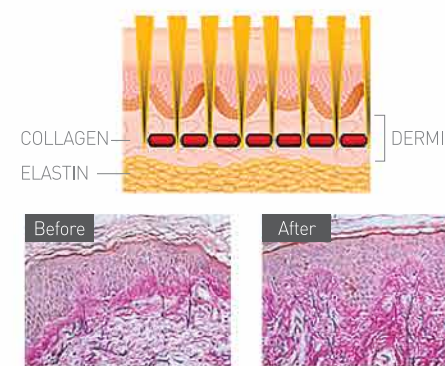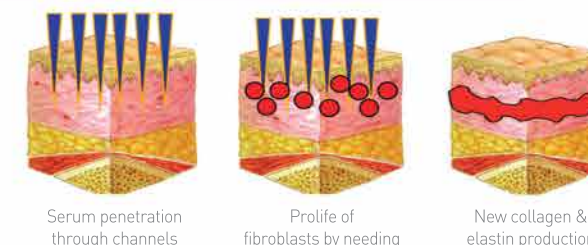

2. WHM(Wound Healing Mechanisms) by RF energy


Since Micro-needles are insulated except the edge part(0.3mm), Micro-needle electrodes of VIVACE deliver bi-polar energy directly deep into dermis without damage in epidermis. Minimally invasive RF thermal energy creates fractional electro coagulation in dermis which facilitate natural Wound Healing Process. WHM promote the remodeling of collagen and wound contraction, which enhance skin laxity. Not only stimulate dermal remodeling but also RF inhibit sebaceous gland activity to improve acne.

Distinction of UIVACE

1. 1MHz/2MHz RF for different target

| | 1MHz | 2MHz |
|------------|---------------------------------------|---|
| Wavelength | | |
| Range | | |
| Indication | Wider RF emission for skin tightening | Intensive RF emission for scar & pore treatment |


2. Provide Insulated & Non-insulated needle

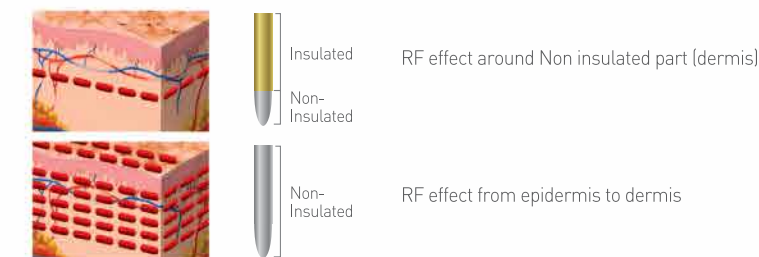

3. Light emitting diode



BLUE LED targets to kill propionibacterium acne or the bacteria associated with acne.

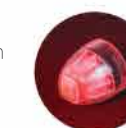

RED LED activates skin cell to increase collagen production moisture levels in skin and improve skin tone & texture.

Clinical Results

